


SALTS ALOE RINGS




MODERATE ABSORPTION AND ADHESION

Features:


- ▶ Skin-friendly, thin and flexible hydrocolloid rings
- ▶ Aloe Vera extracts to help soothe and protect
- ▶ Quick and easy to apply
- ▶ No need for any initial warming or moulding – simply apply straight over the stoma
- ▶ Discreet and low-profile, at only 1.8mm thick
- ▶ Can be stretched to accommodate moderately oval stomas
- ▶ Available in 4 sizes
- ▶ Suitable for all stoma types

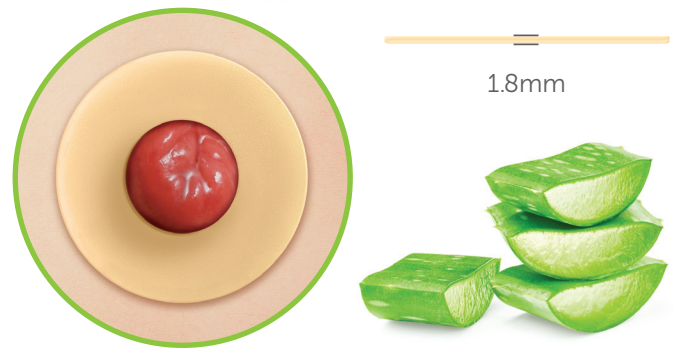


HOW TO APPLY

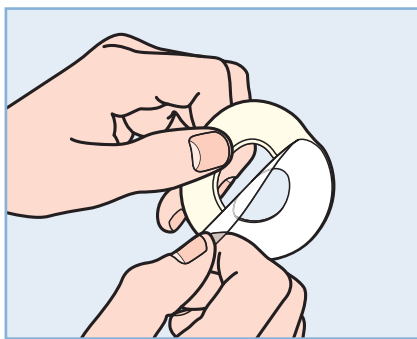

1 Remove the film backing from the ring. No initial warming or moulding is required.

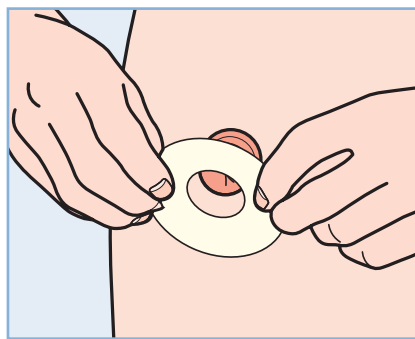

2 Place the ring around the stoma, then smooth round it with the fingers.

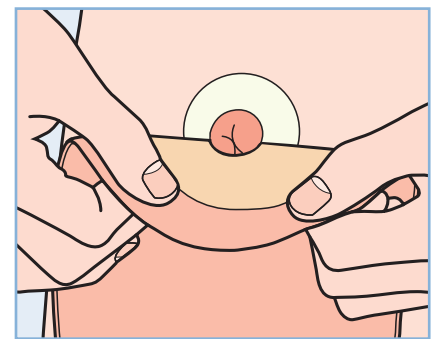

3 Place the pouch gently over the stoma and ring, and smooth all round it with the fingers.

SALTS ALOE RINGS		
INTERNAL DIAMETER	ORDER CODES	PACK SIZE
25mm	SAR25	30
32mm	SAR32	30
38mm	SAR38	30
45mm	SAR45	30